
nanite Documentation

Release 1.5.5

Paul Müller

Jun 05, 2020

1	Getting started	3
1.1	Installation	3
1.2	What is nanite?	3
1.3	Supported file formats	3
1.4	Use cases	4
1.5	Basic usage	4
1.5.1	How to cite	4
2	Command-line interface	5
2.1	nanite-setup-profile	5
2.2	nanite-fit	5
2.3	nanite-rate	5
2.4	nanite-generate-training-set	6
3	Fitting guide	7
3.1	Preprocessors	7
3.2	Models	7
3.3	Parameters	8
3.4	Workflow	8
3.4.1	Command-line usage	8
3.4.2	Scripting usage	10
4	Rating workflow	13
4.1	Rating experimental data manually	13
4.2	Generating the training set	15
4.3	Applying the training set	15
5	Scripting examples	17
5.1	Approximating the Hertzian model with a spherical indenter	17
5.2	Fitting and rating	19
6	Code reference	23
6.1	Module level aliases	23
6.2	Force-indentation data	23
6.3	Groups	26
6.4	Loading data	27
6.5	Preprocessing	27

6.6	Modeling	28
6.6.1	Methods and constants	28
6.6.2	Models	29
6.7	Fitting	35
6.8	Rating	36
6.8.1	Features	36
6.8.2	Rater	39
6.8.3	Regressors	40
6.8.4	Manager	40
6.9	Quantitative maps	42
7	Changelog	45
7.1	version 1.5.5:	45
7.2	version 1.5.4:	45
7.3	version 1.5.3	45
7.4	version 1.5.2	45
7.5	version 1.5.1	46
7.6	version 1.5.0	46
7.7	version 1.4.1	46
7.8	version 1.4.0	46
7.9	version 1.3.0	46
7.10	version 1.2.4	46
7.11	version 1.2.3	46
7.12	version 1.2.2	47
7.13	version 1.2.1	47
7.14	version 1.2.0	47
7.15	version 1.1.2	47
7.16	version 1.1.1	47
7.17	version 1.1.0	47
7.18	version 1.0.1	47
7.19	version 1.0.0	48
7.20	version 0.9.3	48
7.21	version 0.9.2	48
7.22	version 0.9.1	48
7.23	version 0.9.0	48
7.24	version 0.8.0	48
8	Bibliography	49
9	Indices and tables	51
	Bibliography	53
	Python Module Index	55
	Index	57

Nanite is a Python library for loading, fitting, and rating AFM force-distance data of cells and tissues. This is the documentation of nanite version 1.5.5.

1.1 Installation

To install nanite, use one of the following methods (the package dependencies will be installed automatically):

- **from PyPI:** `pip install nanite[CLI]`
- **from sources:** `pip install -e .[CLI]`

The appendix [CLI] makes sure that all dependencies for the *command line interface* are installed. If you are only using nanite as a Python module, you may safely omit it.

Note that if you are installing from source or if no binary wheel is available for your platform and Python version, `Cython` will be installed to build the required nanite extensions. If this process fails, please request a binary wheel for your platform (e.g. Windows 64bit) and Python version (e.g. 3.6) by creating a new [issue](#).

1.2 What is nanite?

The development of nanite was motivated by a unique problem that arises in AFM force-distance data analysis, particularly for biological samples: The data quality varies a lot due to biological variation and due to experimental complexities that have to be dealt with when measuring biological samples. To address this problem, nanite makes use of machine-learning (à la `scikit-learn`), which allows to automatically determine the quality of a force-distance curve based on a user-defined rating scheme (see *Rating workflow* for more information). But nanite is much more than just that. It comes with an extensive set of tools for AFM force-distance data analysis.

1.3 Supported file formats

Nanite relies on the `afmformats` package. A list of supported file formats can be found [here](#).

1.4 Use cases

If you are a frequent AFM user, you might have run into several problems involving data analysis, ranging from simple data fitting to the visualization of quantitative force-distance maps. Here are a few usage examples of nanite:

- You would like to automate your data analysis pipeline from loading force-distance data to displaying a fit to the approach part with a Hertz model for a spherical indenter. You can do so with nanite, either via scripting or via the command-line interface that comes with nanite. For more information, see *Fitting guide*.
- You would like to automatically analyze and visualize maps of force-distance data. This is possible with the `nanite.QMap` class.
- You would like to sort force-distance data according to data quality using your own training set (not the one shipped with nanite). Nanite allows you to create your own training set from your own experimental data, locally. Besides that, you can make use of multiple regressors and visualize the rating e.g. of force-distance maps. For an overview, see *Rating workflow*.

1.5 Basic usage

If you are not interested in scripting, please have a look at the *fitting guide*.

In a Python script, you may use nanite as follows:

```
In [1]: import nanite

In [2]: group = nanite.load_group("data/force-save-example.jpj-force")

In [3]: idnt = group[0] # This group actually as only one indentation curve.

In [4]: idnt.apply_preprocessing(["compute_tip_position",
...:                             "correct_force_offset",
...:                             "correct_tip_offset"])
...:

In [5]: idnt.fit_model(model_key="sneddon_spher")

In [6]: idnt.rate_quality() # 0 means bad, 10 means good quality
Out[6]: 9.060746150910978
```

You can find more examples in the *examples* section.

1.5.1 How to cite

If you use nanite in a scientific publication, please cite Müller et al., *BMC Bioinformatics* (2019) [MAM+19].

Command-line interface

The nanite command-line interface (CLI) simplifies several functionalities of nanite, making fitting, rating, and the generation of training sets accessible to the user.

2.1 nanite-setup-profile

Set up a profile for fitting and rating. The profile is stored in the user's default configuration directory. Setting up a profile is required prior to running *nanite-fit* and *nanite-rate*.

```
usage: nanite-setup-profile [-h]
```

2.2 nanite-fit

Fit AFM force-distance data. Statistics (.tsv file) and visualizations of the fits (multi-page .tif file) are stored in the results directory.

```
usage: nanite-fit [-h] data_path out_dir
```

positional arguments	
data_path	input folder containing AFM force-distance data
out_dir	results directory

2.3 nanite-rate

Manually rate (the fit to) AFM force-distance data. A graphical user interface allows to rate and comment on each force-distance curve. The fits and the raw data are stored in a rating container that can then be passed to *nanite-generate-training-set*.

```
usage: nanite-rate [-h] data_path rating_path
```

positional arguments	
data_path	input folder containing AFM force-distance data
rating_path	path to the output rating container (will be created if it does not already exist)

2.4 nanite-generate-training-set

Create a training set for usage in nanite from rating containers (.h5 files manually created with *nanite-rate*).

```
usage: nanite-generate-training-set [-h] data_path out_dir
```

positional arguments	
data_path	path to a rating container or a folder containing rating containers
out_dir	directory where the training set will be stored

This is a summary of the methods used by nanite for fitting force-distance data. Examples are given below.

3.1 Preprocessors

Prior to data analysis, a force-distance curve has to be preprocessed. One of the most important preprocessing steps is to perform a tip-sample separation which computes the correct tip position from the recorded piezo height and the cantilever deflection. Other preprocessing steps correct for offsets or smoothen the data:

preprocessor key	description	details
compute_tip_position	Compute the tip-sample separation	code reference
correct_force_offset	Correct the force offset with an average baseline value	code reference
correct_split_approach_retract	Split the approach and retract curves (farthest point method)	code reference
correct_tip_offset	Correct the offset of the tip position	code reference
smooth_height	Smoothen height data	code reference

3.2 Models

Nanite comes with a predefined set of model functions that are identified (in scripting as well as in the command line interface) via their model keys.

model key	description	details
hertz_cone	conical indenter (Hertz)	code reference
hertz_para	parabolic indenter (Hertz)	code reference
hertz_pyr3s	pyramidal indenter, three-sided (Hertz)	code reference
sneddon_spher	spherical indenter (Sneddon)	code reference
sneddon_spher_approx	spherical indenter (Sneddon, approximative)	code reference

These model functions can be used to fit experimental force-distance data that have been preprocessed as described above.

3.3 Parameters

Besides the modeling parameters (e.g. Young’s modulus or contact point), nanite allows to define an extensive set of fitting options, that are described in more detail in *nanite.fit.IndentationFitter*.

parameter	description
model_key	Key of the model function used
optimal_fit_edelta	Plateau search for Young’s modulus
optimal_fit_num_samples	Number of points for plateau search
params_initial	Initial parameters
preprocessing	List of preprocessor keys
range_type	‘absolute’ for static range, ‘relative cp’ for dynamic range
range_x	Fitting range (min/max)
segment	Which segment to fit (‘approach’ or ‘retract’)
weight_cp	Suppression of residuals near contact point
x_axis	X-data used for fitting (defaults to ‘top position’)
y_axis	Y-data used for fitting (defaults to ‘force’)

3.4 Workflow

There are two ways to fit force-distance curves with nanite: via the *command line interface (CLI)* or via Python scripting. The CLI does not require programming knowledge while Python-scripting allows fine-tuning and straight-forward automation.

3.4.1 Command-line usage

First, setup up a fitting profile by running (e.g. in a command prompt on Windows).

```
nanite-setup-profile
```

This program will ask you to specify preprocessors, model parameters, and other fitting parameters. Simply enter the values via the keyboard and hit enter to let them be acknowledged. If you want to use the default values, simply hit enter without typing anything. A typical output will look like this:

```
Define preprocessing:
 1: compute_tip_position
 2: correct_force_offset
 3: correct_split_approach_retract
```

(continues on next page)

(continued from previous page)

```

 4: correct_tip_offset
 5: smooth_height
(currently '1,2,4'):

Select model number:
 1: hertz_cone
 2: hertz_para
 3: hertz_pyr3s
 4: sneddon_spher
 5: sneddon_spher_approx
(currently '5'):

Set fit parameters:
- initial value for E [Pa] (currently '3000.0'): 50
  vary E (currently 'True'):
- initial value for R [m] (currently '1e-5'): 18.64e-06
  vary R (currently 'False'):
- initial value for nu (currently '0.5'):
  vary nu (currently 'False'):
- initial value for contact_point [m] (currently '0.0'):
  vary contact_point (currently 'True'):
- initial value for baseline [N] (currently '0.0'):
  vary baseline (currently 'False'):

Select range type (absolute or relative):
(currently 'absolute'):

Select fitting interval:
left [µm] (currently '0.0'):
right [µm] (currently '0.0'):

Suppress residuals near contact point:
size [µm] (currently '0.5'): 2

Select training set:
training set (path or name) (currently 'zef18'):

Select rating regressor:
 1: AdaBoost
 2: Decision Tree
 3: Extra Trees
 4: Gradient Tree Boosting
 5: Random Forest
 6: SVR (RBF kernel)
 7: SVR (linear kernel)
(currently '3'):

Done. You may edit all parameters in '/home/user/.config/nanite/cli_profile.cfg'.

```

In this example, the only modifications of the default values are the initial value of the Young's modulus (50 Pa), the value for the tip radius (18.64 µm), and the suppression of residuals near the contact point with a ± 2 µm interval. When `nanite-setup-profile` is run again, it will use the values from the previous run as default values. The training set and rating regressor options are discussed in the *rating workflow*.

Finally, to perform the actual fitting, use the command-line script

```
nanite-fit data_path output_path
```

This command will recursively search the input folder `data_path` for data files, fit the data with the parameters in the profile, and write the statistics (*statistics.tsv*) and visualizations of the fits (multi-page TIFF file *plots.tif*, open with Fiji or the Windows Photo Viewer) to the directory `output_path`.

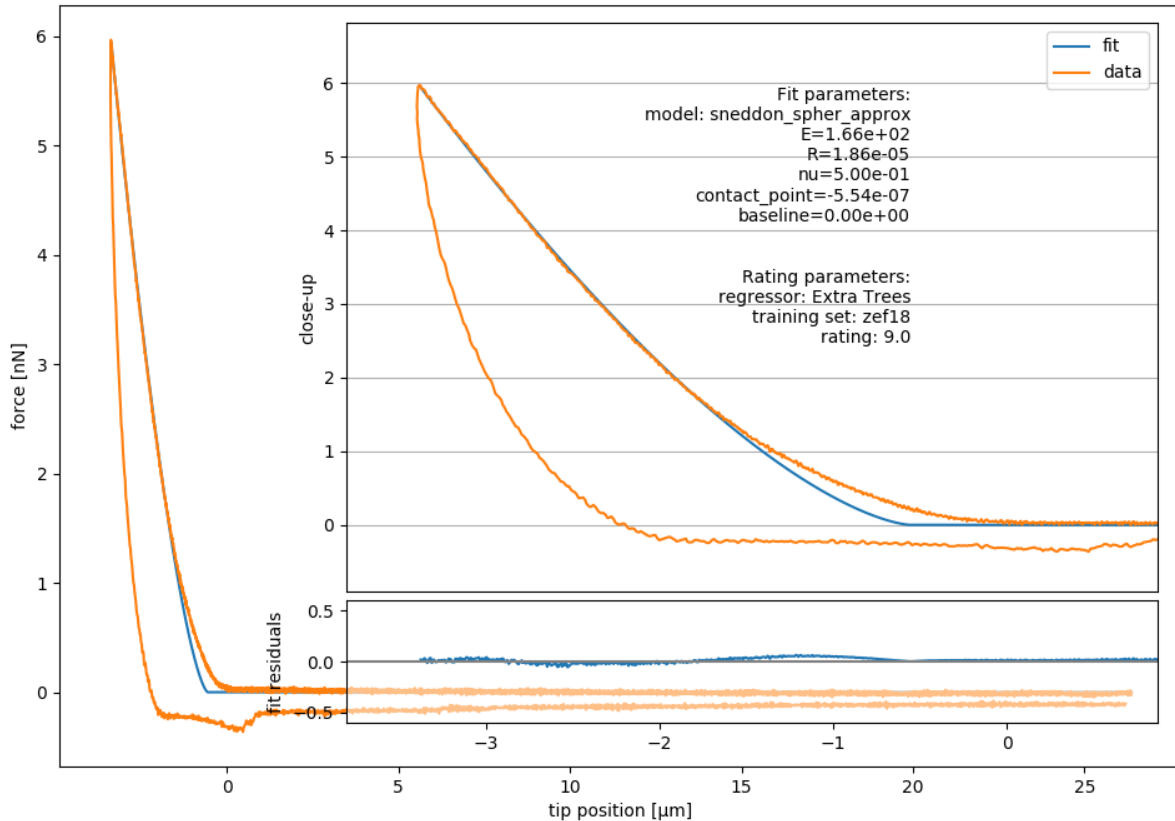


Fig. 3.1: Example image generated with `nanite-fit`. Note that the dataset is already rated with the default method “Extra Trees” and the default training set label “zef18”. See [Rating workflow](#) for more information on rating.

3.4.2 Scripting usage

Using `nanite` in a Python script for data fitting is straight forward. First, load the data; `group` is an instance of `nanite.IndentationGroup`:

```
In [1]: import nanite
In [2]: group = nanite.load_group("data/force-save-example.jpk-force")
```

Second, obtain the first `nanite.Indentation` instance and apply the preprocessing:

```
In [3]: idnt = group[0]
In [4]: idnt.apply_preprocessing(["compute_tip_position",
...:                             "correct_force_offset",
```

(continues on next page)

(continued from previous page)

```
...:         "correct_tip_offset"])
...:
```

Now, setup the model parameters:

```
In [5]: idnt.fit_properties["model_key"] = "sneddon_spher"

In [6]: params = idnt.get_initial_fit_parameters()

In [7]: params["E"].value = 50

In [8]: params["R"].value = 18.64e-06

In [9]: params.pretty_print()
Name          Value      Min      Max  Stderr  Vary  Expr
E              50         0      inf    None   True  None
R          1.864e-05      0      inf    None   False None
baseline        0      -inf    inf    None   True  None
contact_point   0      -inf    inf    None   True  None
nu              0.5         0      0.5    None   False None
```

Finally, fit the model:

```
In [10]: idnt.fit_model(model_key="sneddon_spher", params_initial=params, weight_
↳cp=2e-6)

In [11]: idnt.fit_properties["params_fitted"].pretty_print()
Name          Value      Min      Max  Stderr  Vary  Expr
E            165.8         0      inf  0.1576   True  None
R          1.864e-05      0      inf     0    False  None
baseline   -6.09e-13      -inf    inf  2.316e-13  True  None
contact_point -5.54e-07      -inf    inf  1.416e-09  True  None
nu              0.5         0      0.5     0    False  None
```

The fitting results are identical to those shown in [figure 3.1](#) above.

Note that, amongst other things, preprocessing can also be specified directly in the `fit_model` function.

Rating workflow

One of the main aims of *nanite* is to simplify data analysis by sorting out bad curves automatically based on a user defined rating scheme. *Nanite* allows to automate the rating process using machine learning, based on *scikit-learn*. In short, an estimator is trained with a sample dataset that was manually rated by a user. This estimator is then applied to new data and, in an optimal scenario, reproduces the rating scheme that the user intended when he rated the training dataset. For a more detailed analysis, please refer to [MAM+19].

Nanite already comes with a default training set that is based on AFM data recorded for zebrafish spinal cord sections, called *zef18*. The original *zef18* dataset is available online [MMG18]. Download links:¹

- <https://ndownloader.figshare.com/files/13481393>
- <https://zenodo.org/record/1551200/files/zef18.h5>
- <https://b2share.eudat.eu/api/files/bf481c9b-14ff-47b1-baf5-e569d0199be6/zef18.h5>

With *nanite*, you can also create your own training set. The required steps to do so are described in the following.

4.1 Rating experimental data manually

In the rating step, experimental data are fitted and manually rated by the user. The raw data, the preprocessed data, the fit, all parameters, and the manual rating are then stored in a rating container (an HDF5 file).

First, set up a fitting profile using *nanite-setup-profile* if you have not already done so in the *fitting guide*. You can run the command `nanite-setup-profile` again to verify that all settings are correct.

To start manual rating, use the command *nanite-rate*. The first argument is a folder containing experimental force-distance curves and the second argument is a path to a rating container (`nameXY.h5`). If the rating container already exists, new data will be appended (nothing is overridden).

```
nanite-rate path/to/data/directory path/to/nameXY.h5
```

This will open a graphical user interface that displays the preprocessed and fitted experimental data:

¹ The SHA256 checksum of *zef18.h5* is 63d89a8aa911a255fb4597b2c1801e30ea14810feef1bb42c11ef10f02a1d055.

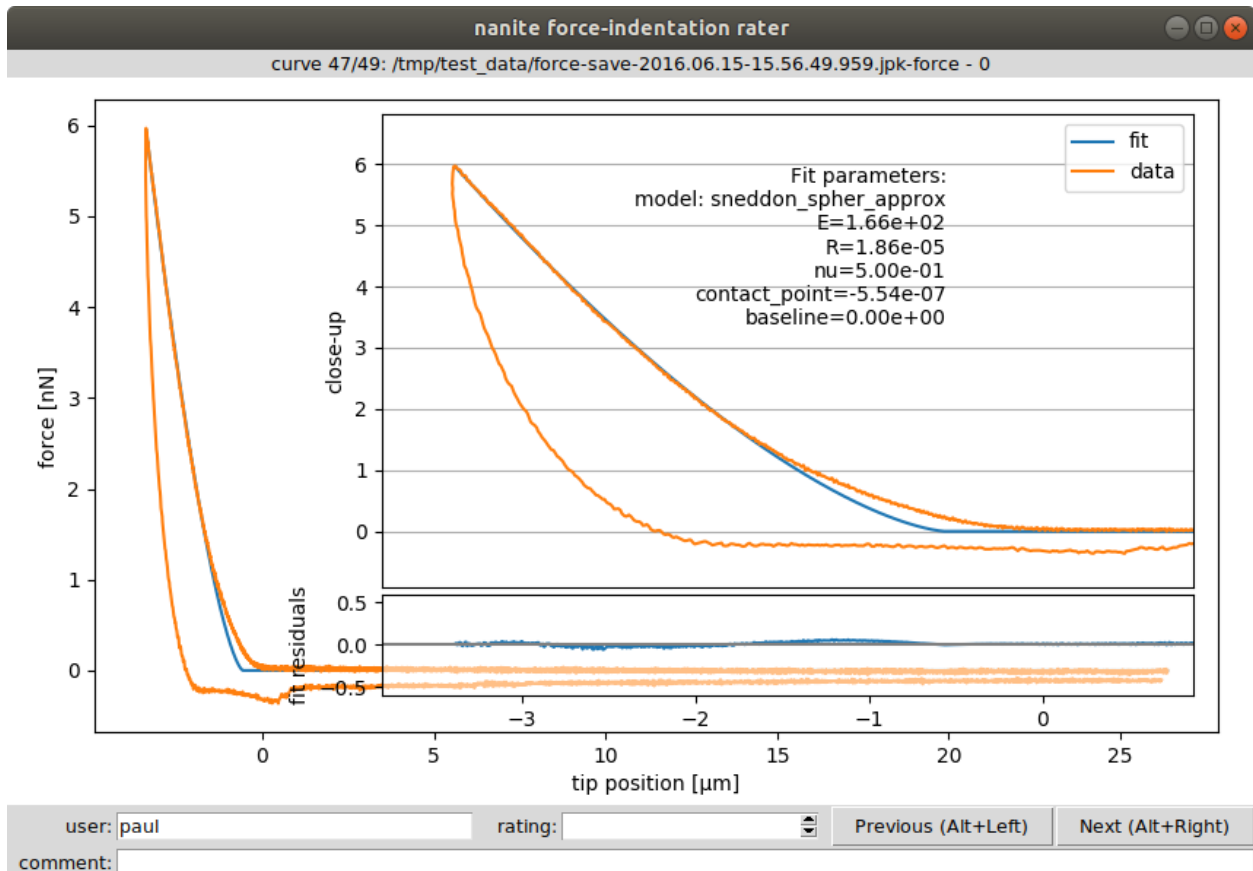


Fig. 4.1: Graphical user interface (GUI) for rating. The inset shows a close-up of the indentation part and the fitted parameters. The user name (defaults to login name) is used to assign a rating to a user (not mandatory). The rating (integer from 0/bad to 10/good or -1/invalid) and a comment can be defined for each curve. The shortcuts `ALT+Left` and `ALT+Right` can be used to navigate within the dataset while keeping the cursor focused in the *rating* field. While navigating, the data are stored in the rating container and the GUI can be closed without data loss.

For the subsequent steps, it is irrelevant whether you create many small rating containers or one global rating container. Many small containers have the advantage that the effect of individual rating sessions could be analyzed separately, while a global rating container keeps all data in one place.

4.2 Generating the training set

The training set consists only of the samples (features of each force-distance curve) and the manual ratings. It is stored as a set of small text files on disk. As described earlier, nanite comes with the predefined *zef18* training set. In this step, a user-defined training set will be generated for use with nanite.

Use the command `nanite-generate-training-set` to convert the rating container(s) to a training set:

```
nanite-generate-training-set path/to/nameXY.h5 path/to/training_set/
```

This will create the folder `path/to/training_set/ts_nameXY` containing several text files, one for each feature and one for the manual rating.

4.3 Applying the training set

To apply the training set when rating curves with `nanite-fit`, you will have to update the profile using `nanite-setup-profile` again (see *fitting guide*). The relevant program output will look like this:

```
[...]
Select training set:
training set (path or name) (currently 'zef18'): path/to/training_set/ts_nameXY

Select rating regressor:
 1: AdaBoost
 2: Decision Tree
 3: Extra Trees
 4: Gradient Tree Boosting
 5: Random Forest
 6: SVR (RBF kernel)
 7: SVR (linear kernel)
(currently '3'):

Done. You may edit all parameters in '/home/user/.config/nanite/cli_profile.cfg'.
```

When running `nanite-fit data_path output_path` now, the new training set is used for rating. The new ratings are stored in `output_path/statistics.tsv` and can be used for further analysis, e.g. quality assessment or sorting.

If you would like to employ a user-defined training set in a Python script, you may do so by specifying the training set path as an argument to `nanite.Indentation.rate_quality`.

5.1 Approximating the Hertzian model with a spherical indenter

There is no closed form for the Hertzian model with a spherical indenter. The force F does not directly depend on the indentation depth δ , but has an indirect dependency via the radius of the circular contact area between indenter and sample a [Sne65]:

$$F = \frac{E}{1 - \nu^2} \left(\frac{R^2 + a^2}{2} \ln \left(\frac{R + a}{R - a} \right) - aR \right)$$

$$\delta = \frac{a}{2} \ln \left(\frac{R + a}{R - a} \right)$$

Here, E is the Young's modulus, R is the radius of the indenter, and ν is the Poisson's ratio of the probed material.

Because of this indirect dependency, fitting this model to experimental data can be time-consuming. Therefore, it is beneficial to approximate this model with a polynomial function around small values of δ/R using the Hertz model for a parabolic indenter as a starting point [Dob18]:

$$F = \frac{4}{3} \frac{E}{1 - \nu^2} \sqrt{R} \delta^{3/2} \left(1 - \frac{1}{10} \frac{\delta}{R} - \frac{1}{840} \left(\frac{\delta}{R} \right)^2 + \frac{11}{15120} \left(\frac{\delta}{R} \right)^3 + \frac{1357}{6652800} \left(\frac{\delta}{R} \right)^4 \right)$$


This example illustrates the error made with this approach. In nanite, the model for a spherical indenter has the identifier `“sneddon_spher”` and the approximate model has the identifier `“sneddon_spher_approx”`.

The plot shows the error for the parabolic indenter model `“hertz_para”` and for the approximation to the spherical indenter model. The maximum indentation depth is set to R . The error made by the approximation of the spherical indenter is more than four magnitudes lower than the maximum force during indentation.

model_spherical_indenter.py

```
1 import matplotlib.pyplot as plt
2 from mpl_toolkits.axes_grid1 import make_axes_locatable
3 from matplotlib.lines import Line2D
4 import matplotlib as mpl
```

(continues on next page)

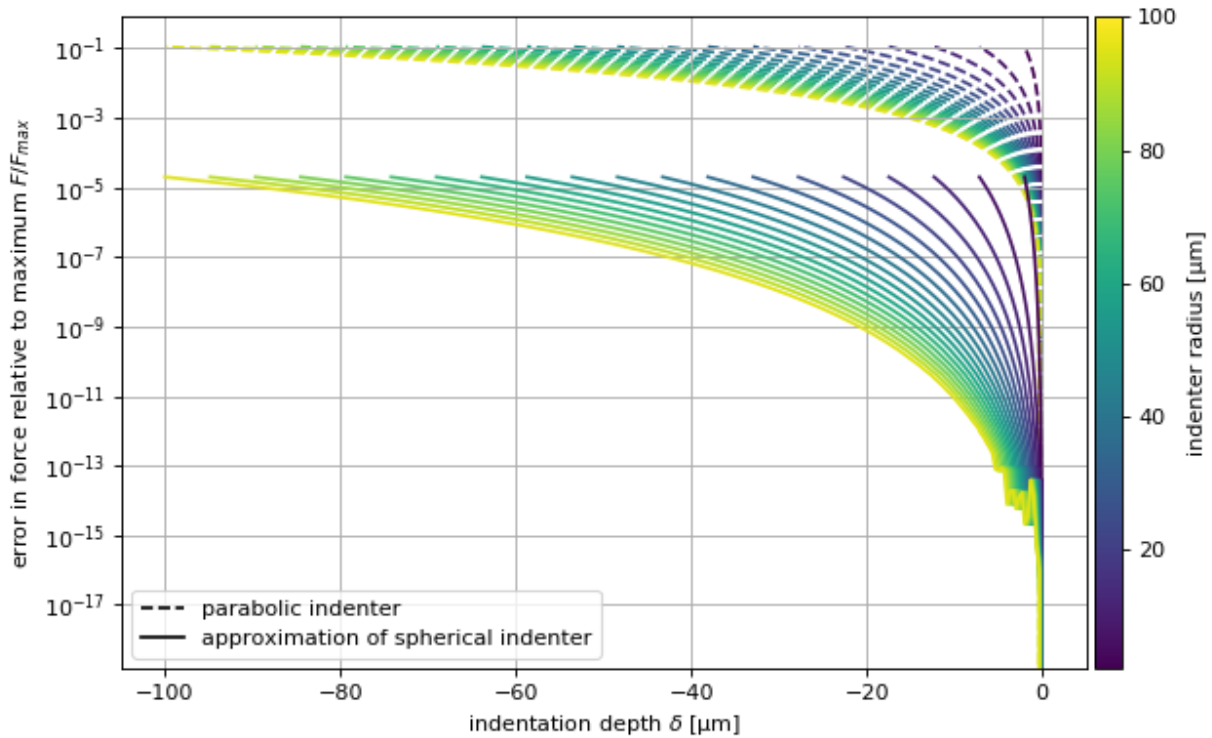

(continued from previous page)

```

5 import numpy as np
6
7 from nanite.model import models_available
8
9 # models
10 exact = models_available["sneddon_spher"]
11 approx = models_available["sneddon_spher_approx"]
12 para = models_available["hertz_para"]
13 # parameters
14 params = exact.get_parameter_defaults()
15 params["E"].value = 1000
16
17 # radii
18 radii = np.linspace(2e-6, 100e-6, 20)
19
20 # plot results
21 plt.figure(figsize=(8, 5))
22
23 # overview plot
24 ax = plt.subplot()
25 for ii, rad in enumerate(radii):
26     params["R"].value = rad
27     # indentation range
28     x = np.linspace(0, -rad, 300)
29     yex = exact.model(params, x)
30     yap = approx.model(params, x)
31     ypa = para.model(params, x)
32     ax.plot(x*1e6, np.abs(yex - yap)/yex.max(),

```

(continues on next page)

(continued from previous page)

```

33         color=mpl.cm.get_cmap("viridis")(ii/radii.size),
34         zorder=2)
35     ax.plot(x*1e6, np.abs(yex - ypa)/yex.max(), ls="--",
36           color=mpl.cm.get_cmap("viridis")(ii/radii.size),
37           zorder=1)
38
39 ax.set_xlabel(r"indentation depth  $\delta$  [ $\mu\text{m}$ "])
40 ax.set_ylabel("error in force relative to maximum  $F/F_{\text{max}}$ ")
41 ax.set_yscale("log")
42 ax.grid()
43
44 # legend
45 custom_lines = [Line2D([0], [0], color="k", ls="--"),
46                Line2D([0], [0], color="k", ls="-"),
47                ]
48 ax.legend(custom_lines, ['parabolic indenter',
49                          'approximation of spherical indenter'])
50
51 divider = make_axes_locatable(ax)
52 cax = divider.append_axes("right", size="3%", pad=0.05)
53
54 norm = mpl.colors.Normalize(vmin=radii[0]*1e6, vmax=radii[-1]*1e6)
55 mpl.colorbar.ColorbarBase(ax=cax,
56                           cmap=mpl.cm.viridis,
57                           norm=norm,
58                           orientation='vertical',
59                           label="indenter radius [ $\mu\text{m}$ "])
60
61
62 plt.tight_layout()
63 plt.show()

```

5.2 Fitting and rating

This example uses a force-distance curve of a zebrafish spinal cord section to illustrate basic data fitting and rating with nanite. The dataset is part of a study on spinal cord stiffness in zebrafish [MKH+19].


fit_and_rate.py

```

1  import matplotlib.gridspec as gridspec
2  import matplotlib.pyplot as plt
3
4  import nanite
5
6  # load the data
7  group = nanite.load_group("data/zebrafish-head-section-gray-matter.jpk-force")
8  idnt = group[0] # this is an instance of `nanite.Indentation`
9  # apply preprocessing
10 idnt.apply_preprocessing(["compute_tip_position",
11                          "correct_force_offset",
12                          "correct_tip_offset"])
13 # set the fit model ("sneddon_spher_approx" is faster than "sneddon_spher"
14 # and sufficiently accurate)
15 idnt.fit_properties["model_key"] = "sneddon_spher_approx"

```

(continues on next page)

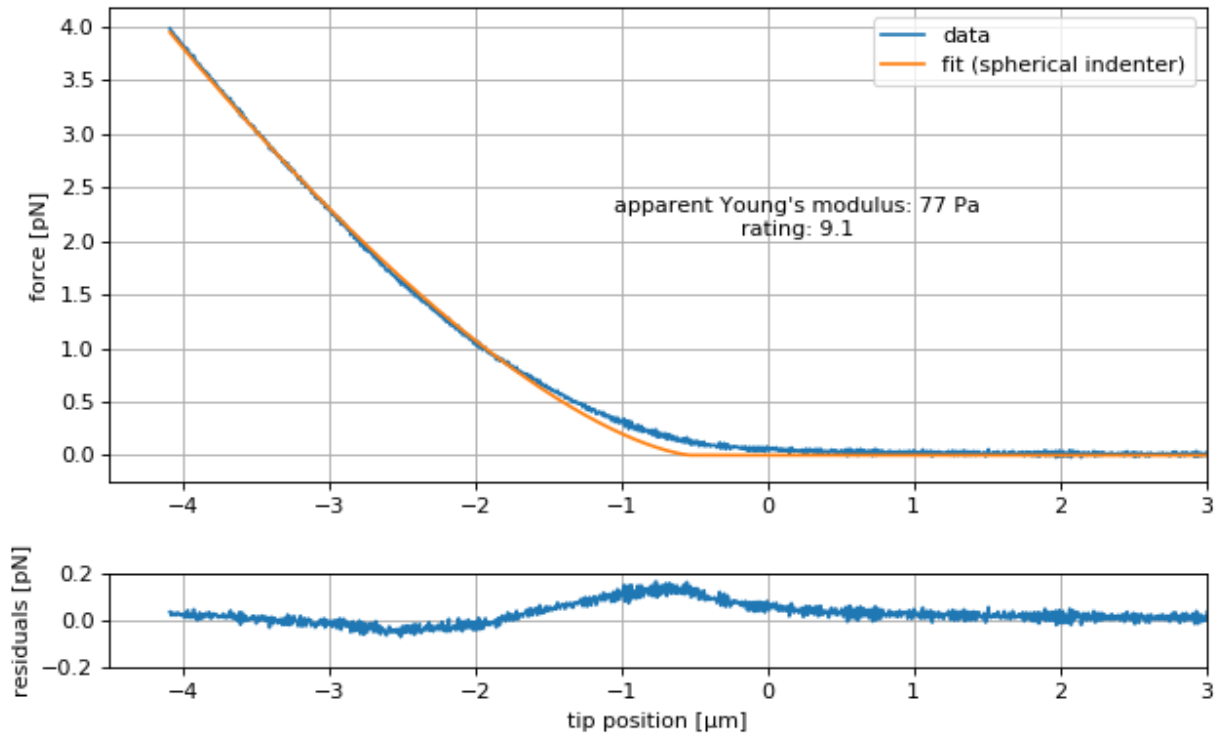

(continued from previous page)

```

16 # get the initial fit parameters
17 params = idnt.get_initial_fit_parameters()
18 # set the correct indenter radius
19 params["R"].value = 18.64e-06
20 # perform the fit with the edited parameters
21 idnt.fit_model(params_initial=params)
22 # obtain the Young's modulus
23 emod = idnt.fit_properties["params_fitted"]["E"].value
24 # obtain a rating for the dataset
25 # (using default regressor and training set)
26 rate = idnt.rate_quality()
27
28 # overview plot
29 plt.figure(figsize=(8, 5))
30 gs = gridspec.GridSpec(2, 1, height_ratios=[5, 1])
31
32 ax1 = plt.subplot(gs[0])
33 ax2 = plt.subplot(gs[1])
34
35 # only plot the approach part (`1` would be retract)
36 where_approach = idnt["segment"] == 0
37
38 # plot force-distance data (nanite uses SI units)
39 ax1.plot(idnt["tip position"][where_approach] * 1e6,
40         idnt["force"][where_approach] * 1e9,
41         label="data")
42 ax1.plot(idnt["tip position"][where_approach] * 1e6,
43         idnt["fit"][where_approach] * 1e9,

```

(continues on next page)

(continued from previous page)

```
44     label="fit (spherical indenter)")
45 ax1.text(.2, 2.05,
46         "apparent Young's modulus: {:.0f} Pa\n".format(emod)
47         + "rating: {:.1f}".format(rate),
48         ha="center")
49 ax1.legend()
50 # plot residuals
51 ax2.plot(idnt["tip position"][where_approach] * 1e6,
52         (idnt["force"] - idnt["fit"])[where_approach] * 1e9)
53
54 # update plot parameters
55 ax1.set_xlim(-4.5, 3)
56 ax1.set_ylabel("force [pN]")
57 ax1.grid()
58 ax2.set_xlim(-4.5, 3)
59 ax2.set_ylim(-.2, .2)
60 ax2.set_ylabel("residuals [pN]")
61 ax2.set_xlabel("tip position [µm]")
62 ax2.grid()
63
64 plt.tight_layout()
65 plt.show()
```


6.1 Module level aliases

For user convenience, the following objects are available at the module level.

```
class nanite.Indentation
    alias of nanite.indent.Indentation

class nanite.IndentationGroup
    alias of nanite.group.IndentationGroup

class nanite.IndentationRater
    alias of nanite.rate.IndentationRater

class nanite.QMap
    alias of nanite.qmap.QMap

nanite.load_group()
    alias of nanite.group.load_group
```

6.2 Force-indentation data

```
class nanite.indent.Indentation (idnt_data)
    Force-indentation

    Parameters idnt_data (nanite.read.IndentationData) – Object holding the experi-
        mental data

apply_preprocessing (preprocessing=None)
    Perform curve preprocessing steps

    Parameters preprocessing (list) – A list of preprocessing method names that are stored
        in the IndentationPreprocessor class. If set to None, self.preprocessing will be used.
```

compute_emodulus_mindelta (*callback=None*)

Elastic modulus in dependency of maximum indentation

The fitting interval is varied such that the maximum indentation depth ranges from the lowest tip position to the estimated contact point. For each interval, the current model is fitted and the elastic modulus is extracted.

Parameters **callback** (*callable*) – A method that is called with the *emoduli* and *indentations* as the computation proceeds every five steps.

Returns **emoduli, indentations** – The fitted elastic moduli at the corresponding maximal indentation depths.

Return type 1d ndarrays

Notes

The information about emodulus and mindelta is also stored in *self.fit_properties* with the keys “optimal_fit_E_array” and “optimal_fit_delta_array”, if *self.fit_model* is called with the argument *search_optimal_fit* set to *True*.

estimate_contact_point_index ()

Estimate the contact point

Contact point (CP) estimation is performed with two methods and that one which returns the smallest index is returned.

Method 1: baseline deviation

1. Obtain the baseline (initial 10% of the approach curve)
2. Compute average and maximum deviation of the baseline
3. The CP is the index of the approach curve where it exceeds twice of the maximum deviation

Method 2: sign of gradient

1. Perform a median filter on the approach curve
2. Compute the gradient
3. Cut off trailing 10 points from the gradient (noise)
4. The CP is the index of the gradient curve when the sign changes, measured from the point of maximal indentation.

If one of the methods fail, the index 0 is returned.

estimate_optimal_mindelta ()

Estimate the optimal indentation depth

This is a convenience function that wraps around *compute_emodulus_mindelta* and *IndentationFitter.compute_opt_mindelta*.

export (*path, fmt='tab'*)

Saves the current data as tab separated values

fit_model (***kwargs*)

Fit the approach-retract data to a model function

Parameters

- **model_key** (*str*) – A key referring to a model in *nanite.model.models_available*

- **params_initial** (*instance of lmfit.Parameters or dict*) – Parameters for fitting. If not given, default parameters are used.
- **range_x** (*tuple of 2*) – The range for fitting, see *range_type* below.
- **range_type** (*str*) – One of:
 - **absolute**: Set the absolute fitting range in values given by the *x_axis*.
 - **relative cp**: In some cases it is desired to be able to fit a model only up until a certain indentation depth (tip position) measured from the contact point. Since the contact point is a fit parameter as well, this requires a two-pass fitting.
- **preprocessing** (*list of str*) – Preprocessing
- **segment** (*str*) – One of “approach” or “retract”.
- **weight_cp** (*float*) – Weight the contact point region which shows artifacts that are difficult to model with e.g. Hertz.
- **optimal_fit_edelta** (*bool*) – Search for the optimal fit by varying the maximal indentation depth and determining a plateau in the resulting Young’s modulus (fitting parameter “E”).

get_ancillary_parameters (*model_key=None*)

Compute ancillary parameters for the current model

get_initial_fit_parameters (*model_key=None, common_ancillaries=True, model_ancillaries=True*)

Return the initial fit parameters

If there are not initial fit parameters set in *self.fit_properties*, then they are computed.

Parameters

- **model_key** (*str*) – Optionally set a model key. This will override the “model_key” key in *self.fit_properties*.
- **global_ancillaries** (*bool*) – Guess global ancillaries such as the contact point.
- **model_ancillaries** (*bool*) – Guess model-related ancillaries

Notes

global_ancillaries and *model_ancillaries* only have an effect if *self.fit_properties*["params_initial"] is set.

get_rating_parameters ()

Return current rating parameters

rate_quality (*regressor='Extra Trees', training_set='zef18', names=None, lda=None*)

Compute the quality of the obtained curve

Uses heuristic approaches to rate a curve.

Parameters

- **regressor** (*str*) – The regressor name used for rating.
- **training_set** (*str*) – A label for a training set shipped with nanite or a path to a training set.

Returns rating – A value between 0 and 10 where 0 is the lowest rating. If no fit has been performed, a rating of -1 is returned.

Return type float

Notes

The rating is cached based on the fitting hash (see *IndentationFitter._hash*).

reset ()

Resets all data operations

data = None

All data as afmformats.AFMForceDistance

fit_properties

Fitting results, see *Indentation.fit_model()*

preprocessing = None

Default preprocessing steps, see *Indentation.apply_preprocessing()*.

6.3 Groups

class nanite.group.IndentationGroup (*path=None, callback=None*)

Group of Indentation

Parameters

- **path** (*str*) – The path to the data file. The data format is determined using the extension of the file and the data is loaded with the correct method.
- **callback** (*callable or None*) – A method that accepts a float between 0 and 1 to externally track the process of loading the data.

append (*item*)

Append an indentation dataset

Parameters **item** (*nanite.indent.Indentation*) – Force-indentation dataset

get_enum (*enum*)

Return the indentation curve with this enum value

Raises

- ValueError if multiple curves with the same enum value exist.
- KeyError if the enum value is not found

index (*item*)

subgroup_with_path (*path*)

Return a subgroup with measurements matching *path*

nanite.group.load_group (*path, callback=None*)

Load indentation data from disk

Parameters

- **path** (*path-like*) – Path to experimental data
- **callback** (*callable or None*) – Callback function for tracking loading progress

Returns **group** – Indentation group with force-distance data

Return type nanite.IndentationGroup

6.4 Loading data

`nanite.read.get_data_paths` (*path*)

Obtain a list of data files

Parameters *path* (*str* or *pathlib.Path*) – Path to a data file or a directory containing data files.

Returns *paths* – All supported data files found in *path*. If *path* is a file, [*pathlib.Path(path)*] is returned. If *path* has an unsupported extion, an empty list is returned.

Return type list of *pathlib.Path*

`nanite.read.get_data_paths_enum` (*path*, *skip_errors=False*)

`nanite.read.load_data` (*path*, *callback=None*)

Load data and return list of *afmformats.AFMForceDistance*

6.5 Preprocessing

exception `nanite.preproc.CannotSplitWarning`

class `nanite.preproc.IndentationPreprocessor`

static apply (*apret*, *preproc_names*)

Perform force-distance preprocessing steps

Parameters

- **apret** (`nanite.Indentation`) – The afm data to preprocess
- **preproc_names** (*list*) – A list of names for static methods in *IndentationPreprocessor* that will be applied (in the order given).

Notes

This method is usually called from within the *Indentation* class instance. If you are using this class directly and apply it more than once, you might need to call *apret.reset()* before preprocessing a second time.

static available ()

List available preprocessor names

static compute_tip_position (*apret*)

Compute the tip-sample separation

This computation correctly reproduces the column “Vertical Tip Position” as it is exported by the JPK analysis software with the checked option “Use Unsmoothed Height”.

static correct_force_offset (*apret*)

Correct the force offset with an average baseline value

static correct_split_approach_retract (*apret*)

Split the approach and retract curves (farthest point method)

Approach and retract curves are defined by the microscope. When the direction of piezo movement is flipped, the force at the sample tip is still increasing. This can be either due to a time lag in the AFM system or due to a residual force acting on the sample due to the bent cantilever.

To repair this time lag, we append parts of the retract curve to the approach curve, such that the curves are split at the minimum height.

static correct_tip_offset (*apret*)

Correct the offset of the tip position

An estimate of the tip position is used to compute the contact point.

static smooth_height (*apret*)

Smoothen height data

For the columns “height (measured)” and “tip position”, and for the approach and retract data separately, this method adds the columns “height (measured, smoothed)” and “tip position (smoothed)” to *self.data*.

```
nanite.preproc.available_preprocessors = ['compute_tip_position', 'correct_force_offset',  
Available preprocessors
```

6.6 Modeling

6.6.1 Methods and constants

exception `nanite.model.ModelIncompleteError`

`nanite.model.get_anc_parm_keys` (*model_key*)

Return the key names of a model’s ancillary parameters

`nanite.model.get_anc_parms` (*idnt, model_key*)

Compute ancillary parameters for a force-distance dataset

Ancillary parameters include parameters that:

- are unrelated to fitting: They may just be important parameters to the user.
- require the entire dataset: They cannot be extracted during fitting, because they require more than just the approach xor retract curve to compute (e.g. hysteresis, jump of retract curve at maximum indentation). They may, additionally, depend on initial fit parameters set by the user.
- require a fit: They are dependent on fitting parameters but are not required during fitting.

Notes

If an ancillary parameter name matches that of a fitting parameter, then it is assumed that it can be used for fitting. Please see `nanite.indent.Indentation.get_initial_fit_parameters()` and `nanite.fit.guess_initial_parameters()`.

Ancillary parameters are set to *np.nan* if they cannot be computed.

Parameters

- **idnt** (`nanite.indent.Indentation`) – The force-distance data for which to compute the ancillary parameters
- **model_key** (*str*) – Name of the model

Returns **ancillaries** – key-value dictionary of ancillary parameters

Return type `collections.OrderedDict`

`nanite.model.get_init_parms` (*model_key*)

Get initial fit parameters for a model

`nanite.model.get_model_by_name(name)`

Convenience function to obtain a model by name instead of by key

`nanite.model.get_parm_name(model_key, parm_key)`

Return parameter label

Parameters

- **model_key** (*str*) – The model key (e.g. “hertz_cone”)
- **parm_key** (*str*) – The parameter key (e.g. “E”)

Returns **parm_name** – The parameter label (e.g. “Young’s Modulus”)

Return type *str*

`nanite.model.get_parm_unit(model_key, parm_key)`

Return parameter unit

Parameters

- **model_key** (*str*) – The model key (e.g. “hertz_cone”)
- **parm_key** (*str*) – The parameter key (e.g. “E”)

Returns **parm_unit** – The parameter unit (e.g. “Pa”)

Return type *str*

`nanite.model.register_model(module, module_name)`

Register a fitting model

`nanite.model.ANCILLARY_COMMON = {'max_indent': ('Maximum indentation', 'm')}`

Common ancillary parameters

6.6.2 Models

Each model is implemented as a submodule in *nanite.model*. For instance `nanite.model.model_hertz_parabolic`. Each of these modules implements the following functions (which are not listed for each model in the subsections below):

`nanite.model.model_submodule.get_parameter_defaults()`

Return the default parameters of the model.

`nanite.model.model_submodule.model()`

Wrap the actual model for fitting.

`nanite.model.model_submodule.residual()`

Compute the residuals during fitting.

In addition, each submodule contains the following attributes:

`nanite.model.model_submodule.model_doc`

The doc-string of the model function.

`nanite.model.model_submodule.model_key`

The model key used in the command line interface and during scripting.

`nanite.model.model_submodule.model_name`

The name of the model.

`nanite.model.model_submodule.parameter_keys`

Parameter keywords of the model for higher-level applications.

`nanite.model.model_submodule.parameter_names`

Parameter names of the model for higher-level applications.

`nanite.model.model_submodule.parameter_units`

Parameter units for higher-level applications.

conical indenter (Hertz)

model key	hertz_cone
model name	conical indenter (Hertz)
model location	nanite.model.model_conical_indenter

`nanite.model.model_conical_indenter.hertz_conical`(*E*, *delta*, *alpha*, *nu*, *contact_point=0*, *baseline=0*)

Hertz model for a conical indenter

$$F = \frac{2 \tan \alpha}{\pi} \frac{E}{1 - \nu^2} \delta^2$$

Parameters

- **E** (*float*) – Young’s modulus [N/m²]
- **delta** (*1d ndarray*) – Indentation [m]
- **alpha** (*float*) – Half cone angle [degrees]
- **nu** (*float*) – Poisson’s ratio
- **contact_point** (*float*) – Indentation offset [m]
- **baseline** (*float*) – Force offset [N]
- **negindent** (*bool*) – If *True*, will assume that the indentation value(s) given by *delta* are negative and must be multiplied by -1.

Returns **F** – Force [N]

Return type `float`

Notes

These approximations are made by the Hertz model:

- The sample is isotropic.
- The sample is a linear elastic solid.
- The sample is extended infinitely in one half space.
- The indenter is not deformable.
- There are no additional interactions between sample and indenter.

Additional assumptions:

- infinitely sharp probe

References

Love (1939) [Lov39]

parabolic indenter (Hertz)

model key	hertz_para
model name	parabolic indenter (Hertz)
model location	nanite.model.model_hertz_paraboloidal

`nanite.model.model_hertz_paraboloidal.hertz_paraboloidal` (E , δ , R , ν , $\text{contact_point}=0$, $\text{baseline}=0$)

Hertz model for a paraboloidal indenter

$$F = \frac{4}{3} \frac{E}{1 - \nu^2} \sqrt{R} \delta^{3/2}$$

Parameters

- **E** (*float*) – Young’s modulus [N/m²]
- **delta** (*1d ndarray*) – Indentation [m]
- **R** (*float*) – Tip radius [m]
- **nu** (*float*) – Poisson’s ratio
- **contact_point** (*float*) – Indentation offset [m]
- **baseline** (*float*) – Force offset [N]
- **negindent** (*bool*) – If *True*, will assume that the indentation value(s) given by *delta* are negative and must be multiplied by -1.

Returns **F** – Force [N]

Return type `float`

Notes

The original model reads

$$F = \frac{4}{3} \frac{E}{1 - \nu^2} \sqrt{2k} \delta^{3/2},$$

where k is defined by the paraboloid equation

$$\rho^2 = 4kz.$$

These approximations are made by the Hertz model:

- The sample is isotropic.
- The sample is a linear elastic solid.
- The sample is extended infinitely in one half space.
- The indenter is not deformable.
- There are no additional interactions between sample and indenter.

Additional assumptions:

- radius of spherical cell is larger than the indentation

References

Sneddon (1965) [Sne65]

pyramidal indenter, three-sided (Hertz)

model key	hertz_pyr3s
model name	pyramidal indenter, three-sided (Hertz)
model location	nanite.model.model_hertz_three_sided_pyramid

nanite.model.model_hertz_three_sided_pyramid.**hertz_three_sided_pyramid**(*E*,
delta,
al-
pha,
nu,
con-
tact_point=0,
base-
line=0)

Hertz model for three sided pyramidal indenter

$$F = 0.887 \tan \alpha \cdot \frac{E}{1 - \nu^2} \delta^2$$

Parameters

- **E** (*float*) – Young’s modulus [N/m²]
- **delta** (*1d ndarray*) – Indentation [m]
- **alpha** (*float*) – Inclination angle of the pyramidal face [degrees]
- **nu** (*float*) – Poisson’s ratio
- **contact_point** (*float*) – Indentation offset [m]
- **baseline** (*float*) – Force offset [N]
- **negindent** (*bool*) – If *True*, will assume that the indentation value(s) given by *delta* are negative and must be multiplied by -1.

Returns **F** – Force [N]

Return type *float*

Notes

These approximations are made by the Hertz model:

- The sample is isotropic.
- The sample is a linear elastic solid.
- The sample is extended infinitely in one half space.
- The indenter is not deformable.
- There are no additional interactions between sample and indenter.
- The inclination angle of the pyramidal face (in radians) must be close to zero.

References

Bilodeau et al. 1992 [Bil92]

spherical indenter (Sneddon)

model key	sneddon_spher
model name	spherical indenter (Sneddon)
model location	nanite.model.model_sneddon_spherical

`nanite.model.model_sneddon_spherical.delta_of_a()`

Compute indentation from contact area radius (wrapper)

`nanite.model.model_sneddon_spherical.get_a()`

Compute the contact area radius (wrapper)

`nanite.model.model_sneddon_spherical.hertz_spherical()`

Hertz model for Spherical indenter - modified by Sneddon

$$F = \frac{E}{1 - \nu^2} \left(\frac{R^2 + a^2}{2} \ln \left(\frac{R + a}{R - a} \right) - aR \right)$$

$$\delta = \frac{a}{2} \ln \left(\frac{R + a}{R - a} \right)$$

(a is the radius of the circular contact area between bead and sample.)

Parameters

- **E** (*float*) – Young’s modulus [N/m²]
- **delta** (*1d ndarray*) – Indentation [m]
- **R** (*float*) – Tip radius [m]
- **nu** (*float*) – Poisson’s ratio
- **contact_point** (*float*) – Indentation offset [m]
- **baseline** (*float*) – Force offset [N]
- **negindent** (*bool*) – If *True*, will assume that the indentation value(s) given by *delta* are negative and must be multiplied by -1.

Returns **F** – Force [N]

Return type `float`

Notes

These approximations are made by the Hertz model:

- The sample is isotropic.
- The sample is a linear elastic solid.
- The sample is extended infinitely in one half space.
- The indenter is not deformable.
- There are no additional interactions between sample and indenter.

Additional assumptions:

- no surface forces

References

Sneddon (1965) [Sne65]

spherical indenter (Sneddon, approximative)

model key	sneddon_spher_approx
model name	spherical indenter (Sneddon, approximative)
model location	nanite.model.model_sneddon_spherical_approximation

nanite.model.model_sneddon_spherical_approximation.**hertz_sneddon_spherical_approx**(*E*,
delta,
R,
nu,
contact_point=0,
baseline=0)

Hertz model for Spherical indenter - approximation

$$F = \frac{4}{3} \frac{E}{1 - \nu^2} \sqrt{R} \delta^{3/2} \left(1 - \frac{1}{10} \frac{\delta}{R} - \frac{1}{840} \left(\frac{\delta}{R} \right)^2 + \frac{11}{15120} \left(\frac{\delta}{R} \right)^3 + \frac{1357}{6652800} \left(\frac{\delta}{R} \right)^4 \right)$$

Parameters

- **E** (*float*) – Young’s modulus [N/m²]
- **delta** (*1d ndarray*) – Indentation [m]
- **R** (*float*) – Tip radius [m]
- **nu** (*float*) – Poisson’s ratio
- **contact_point** (*float*) – Indentation offset [m]
- **baseline** (*float*) – Force offset [N]
- **negindent** (*bool*) – If *True*, will assume that the indentation value(s) given by *delta* are negative and must be multiplied by -1.

Returns **F** – Force [N]

Return type *float*

Notes

These approximations are made by the Hertz model:

- The sample is isotropic.
- The sample is a linear elastic solid.
- The sample is extended infinitely in one half space.

- The indenter is not deformable.
- There are no additional interactions between sample and indenter.

Additional assumptions:

- no surface forces

References

Sneddon (1965) [Sne65], Dobler (personal communication, 2018) [Dob18]

6.7 Fitting

exception `nanite.fit.FitDataError`

exception `nanite.fit.FitKeyError`

exception `nanite.fit.FitWarning`

class `nanite.fit.FitProperties`

Fit property manager class

Provide convenient access to fit properties as a dictionary and dynamically manage resets due to new initial parameters.

Dynamic properties include:

- set “`params_initial`” to `None` if the “`model_key`” changes
- remove all keys except those in `FP_DEFAULT` if a key that is in `FP_DEFAULT` changes (All other keys are considered to be obsolete fitting results).

Additional attributes:

- “`segment_bool`”: **bool** `False` for “approach” and `True` for “retract”

reset ()

restore (*props*)

update the dictionary without removing any keys

class `nanite.fit.IndentationFitter` (*idnt*, ***kwargs*)

Fit force-distance curves

Parameters

- **idnt** (`nanite.indent.Indentation`) – The dataset to fit
- **model_key** (*str*) – A key referring to a model in `nanite.model.models_available`
- **params_initial** (*instance of `lmfit.Parameters`*) – Parameters for fitting. If not given, default parameters are used.
- **range_x** (*tuple of 2*) – The range for fitting, see *range_type* below.
- **range_type** (*str*) – One of:
 - **absolute**: Set the absolute fitting range in values given by the *x_axis*.
 - **relative cp**: In some cases it is desired to be able to fit a model only up until a certain indentation depth (tip position) measured from the contact point. Since the contact point is a fit parameter as well, this requires a two-pass fitting.

- **preprocessing** (*list of str*) – Preprocessing
- **segment** (*str*) – One of “approach” or “retract”.
- **weight_cp** (*float*) – Weight the contact point region which shows artifacts that are difficult to model with e.g. Hertz.
- **optimal_fit_edelta** (*bool*) – Search for the optimal fit by varying the maximal indentation depth and determining a plateau in the resulting Young’s modulus (fitting parameter “E”).
- **optimal_fit_num_samples** (*int*) – Number of samples to use for searching the optimal fit

compute_emodulus_vs_mindelta (*callback=None*)
Compute elastic modulus vs. minimal indentation curve

static compute_opt_mindelta (*emoduli, indentations*)
Determine the plateau of an emodulus-indentation curve

The following procedure is performed:

1. Smooth the emodulus data with a Butterworth filter
2. Label sequences that have similar values by binning into ten regions between the min and max.
3. Ignore sequences with emodulus that is smaller than the binning size.
4. Determine the longest sequence.

fit ()
Fit the approach-retract data to a model function

get_initial_parameters (*idnt=None, model_key='hertz_para'*)
Get initial fit parameters for a specific model

`nanite.fit.guess_initial_parameters` (*idnt=None, model_key='hertz_para', common_ancillaries=True, model_ancillaries=True*)
Guess initial fitting parameters

Parameters

- **idnt** (`nanite.indent.Indentation`) – The dataset to use for guessing initial fitting parameters using ancillary parameters
- **model_key** (*str*) – The model key
- **common_ancillaries** (*bool*) – Guess global ancillary parameters (such as contact point)
- **model_ancillaries** (*bool*) – Use model-related ancillary parameters

`nanite.fit.obj2str` (*obj*)
String representation of an object for hashing

6.8 Rating

6.8.1 Features

class `nanite.rate.features.IndentationFeatures` (*dataset=None*)

static compute_features (*idnt*, *which_type='all'*, *names=None*, *ret_names=False*)

Compute the features for a data set

Parameters

- **idnt** (*nanite.Indentation*) – A dataset to rate
- **names** (*list of str*) – The names of the rating methods to use, e.g. ["rate_apr_bumps", "rate_apr_mon_incr"]. If None (default), all available rating methods are used.

Notes

names may include features that are excluded by *which_type*. E.g. if a “bool” feature is in *names* but *which_type* is “float”, then the “bool” feature will be silently ignored.

feat_bin_apr_spikes_count ()

spikes during IDT

Sudden spikes in indentation curve

feat_bin_cp_position ()

CP outside of data range

Contact point position outside of range

feat_bin_size ()

dataset too small

Number of points in indentation curve

feat_con_apr_flatness ()

flatness of APR residuals

fraction of the positive-gradient residuals in the approach part

feat_con_apr_size ()

relative APR size

length of the approach part relative to the indentation part

feat_con_apr_sum ()

residuals of APR

absolute sum of the residuals in the approach part

feat_con_bln_slope ()

slope of BLN

slope obtained from a linear least-squares fit to the baseline

feat_con_bln_variation ()

variation in BLN

comparison of the forces at the beginning and at the end of the baseline

feat_con_cp_curvature ()

curvature at CP

curvature of the force-distance data at the contact point

feat_con_cp_magnitude ()

residuals at CP

mean value of the residuals around the contact point

feat_con_idt_maxima_75perc()
maxima in IDT residuals

sum of the indentation residuals' maxima in three intervals in-between 25% and 100% relative to the maximum indentation

feat_con_idt_monotony()
monotony of IDT

change of the gradient in the indentation part

feat_con_idt_spike_area()
area of IDT spikes

area of spikes appearing in the indentation part

feat_con_idt_sum()
overall IDT residuals

sum of the residuals in the indentation part

feat_con_idt_sum_75perc()
residuals in 75% IDT

sum of the residuals in the indentation part in-between 25% and 100% relative to the maximum indentation

classmethod get_feature_funcs (*which_type='all', names=None*)
Return functions that compute features from a dataset

Parameters

- **names** (*list of str*) – The names of the rating methods to use, e.g. ["rate_apr_bumps", "rate_apr_mon_incr"]. If None (default), all available rating methods are returned.
- **which_type** (*str*) – Which features to return: ["all", "bool", "float"].

Returns *raters* – Each item in the list consists contains the name of the rating method and the corresponding rating method.

Return type list of tuples (name, callable)

classmethod get_feature_names (*which_type='all', names=None, ret_indices=False*)
Get features names

Parameters

- **which_type** (*str or list of str*) – Return only features that are of a certain type. See *VALID_FEATURE_TYPES* for valid strings.
- **names** (*list of str*) – Return only features that are in this list.
- **ret_indices** (*bool*) – If True, also return the internal feature indices.

Returns *name_list* – List of feature names (callables of this class)

Return type list of str

contact_point

datafit_apr

datares_apr

dataset = None

current dataset from which features are computed

datax_apr
datay_apr
has_contact_point
is_fitted
is_valid
meta

`nanite.rate.features.VALID_FEATURE_TYPES = ['all', 'binary', 'continuous']`
 Valid keyword arguments for feature types

6.8.2 Rater

class `nanite.rate.rater.IndentationRater` (*regressor=None, scale=None, lda=None, training_set=None, names=None, weight=True, sample_weight=None, *args, **kwargs*)

Rate quality

Parameters

- **regressor** (*sciki-learn RegressorMixin*) – The regressor used for rating
- **scale** (*bool*) – If True, apply a Standard Scaler. If a regressor based on decision trees is used, the Standard Scaler is not used by default, otherwise it is.
- **lda** (*bool*) – If True, apply a Linear Discriminant Analysis (LDA). If a regressor based on a decision tree is used, LDA is not used by default, otherwise it is.
- **training_set** (*tuple of (X, y)*) – The training set (samples, response)
- **names** (*list of str*) – Feature names to use
- **weight** (*bool*) – Weight the input samples by the number of occurrences or with *sample_weight*. For tree-based classifiers, set this to True to avoid bias.
- **sample_weight** (*list-like*) – The sample weights. If set to *None* sample weights are computed from the training set.
- ***args** (*list*) – Positional arguments for `IndentationFeatures`
- ****kwargs** – Keyword arguments for `IndentationFeatures`

See also:

`sklearn.preprocessing.StandardScaler` Standard scaler

`sklearn.discriminant_analysis.LinearDiscriminantAnalysis` Linear discriminant analysis

`nanite.rate.regressors.reg_trees` List of regressors that are identified as tree-based

static `compute_sample_weight` (*X, y*)
 Weight samples according to occurrence in *y*

static `get_training_set_path` (*label='zef18'*)
 Return the path to a training set shipped with nanite

Training sets are stored in the `nanite.rate` module path with `ts_` prepended to *label*.

classmethod `load_training_set` (*path=None, names=None, which_type=['continuous'], remove_nan=True, ret_names=False*)

Load a training set from a directory

By default, only the “continuous” features are imported. The “binary” features are not needed for training; they are used to sort out new force-distance data.

rate (*samples=None, datasets=None*)

Perform rating step

Parameters

- **samples** (*1d or 2d ndarray (cast to 2d ndarray) or None*) – Measured samples, if set to *None*, *dataset* must be given.
- **dataset** (*list of nanite.Indentation*) – Full, fitted measurement

Returns `ratings` – Resulting ratings

Return type `list`

names = None

feature names used by the regressor pipeline

pipeline = None

sklearn pipeline with transforms (and regressor if given)

`nanite.rate.rater.get_available_training_sets()`

List of internal training sets

`nanite.rate.rater.get_rater` (*regressor, training_set='zef18', names=None, lda=None, **reg_kwargs*)

Convenience method to get a rater

Parameters

- **regressor** (*str or RegressorMixin*) – If a string, must be in *reg_names*.
- **training_set** (*str or pathlib.Path or tuple (X, y)*) – A string label representing a training set shipped with nanite, the path to a training set, or a tuple representing the training set (samples, response) for use with sklearn.

Returns `irater` – The rating instance.

Return type `nanite.IndentationRater`

6.8.3 Regressors

scikit-learn regressors and their keyword arguments

`nanite.rate.regressors.reg_names = ['AdaBoost', 'Decision Tree', 'Extra Trees', 'Gradient Boosting Machine']`
List of available default regressor names

`nanite.rate.regressors.reg_trees = ['AdaBoostRegressor', 'DecisionTreeRegressor', 'ExtraTreeRegressor']`
List of tree-based regressor class names (used for keyword defaults in `IndentationRater`)

6.8.4 Manager

Save and load user-rated datasets

class `nanite.rate.io.RateManager` (*path, verbose=0*)

Manage user-defined rates

`export_training_set` (*path*)

`get_cross_validation_score` (*regressor, training_set=None, n_splits=20, random_state=42*)
Regressor cross-validation scoring

Cross-validation is used to identify regressors that over-fit the train set by splitting the train set into multiple learn/test sets and quantifying the regressor performance for each split.

Parameters

- **regressor** (*str* or *RegressorMixin*) – If a string, must be in *reg_names*.
- **training_set** (*X, y*) – If given, do not use `self.samples`

Notes

A `skimage.model_selection.KFold` cross validator is used in combination with the mean squared error score.

Cross-validation score is computed from samples that are filtered with the binary features and only from samples that do not contain any nan values.

`get_rates` (*which='user', training_set='zef18'*)

which: *str* Which rating to return: “user” or a regressor name

`get_training_set` (*which_type='all', prefilter_binary=False, remove_nans=False, transform=False*)
Return (*X, y*) training set

datasets

path = None

Path to the manual ratings (directory or .h5 file)

ratings

samples

The individual sample ratings computed by `afmlib`

verbose = None

verbosity level

`nanite.rate.io.hdf5_rated` (*h5path, indent*)

Test whether an indentation has already been rated

Returns

Return type `is_rated, rating, comment`

`nanite.rate.io.load` (*path, meta_only=False, verbose=0*)

Notes

The `.fit_properties` attribute of each `Indentation` instance is overridden by a simple dictionary, so its functionalities are not available anymore.

`nanite.rate.io.load_hdf5` (*path, meta_only=False*)

`nanite.rate.io.save_hdf5` (*h5path, indent, user_rate, user_name, user_comment, h5mode='a'*)

Store all relevant data of a user rating into an hdf5 file

Parameters

- **h5path** (*str*) – Path to HDF5 rating container where data will be stored
- **indent** (`nanite.Indentation`) – The experimental data processed and fitted with nanite
- **user_rate** (*float*) – Rating given by the user
- **user_name** (*str*) – Name of the rating user

`nanite.rate.io.hash_file`

Compute sha256 hex-hash of a file

Parameters

- **path** (*str or pathlib.Path*) – path to the file
- **blocksize** (*int*) – block size read from the file

Returns `hex` – The first six characters of the hash

Return type `str`

6.9 Quantitative maps

exception `nanite.qmap.DataMissingWarning`

class `nanite.qmap.QMap` (*path_or_dataset, callback=None*)

Quantitative force spectroscopy map handling

Parameters

- **path_or_dataset** (*str or nanite.IndentationGroup*) – The path to the data file. The data format is determined using the extension of the file and the data is loaded with the correct method.
- **callback** (*callable or None*) – A method that accepts a float between 0 and 1 to externally track the process of loading the data.

`feat_data_min_height_measured_um` (*idnt*)

`feat_fit_contact_point` (*idnt*)

`feat_fit_youngs_modulus` (*idnt*)

`feat_meta_rating` (*idnt*)

`feat_meta_scan_order` (*idnt*)

`get_qmap` (*feature, qmap_only=False*)

Return the quantitative map for a feature

Parameters

- **feature** (*str*) – Feature to compute map for (see `QMap.features`)
- **qmap_only** – Only return the quantitative map data, not the coordinates

Returns

- **x, y** (*1d ndarray*) – Only returned if `qmap_only` is False; Pixel grid coordinates along x and y
- **qmap** (*2d ndarray*) – Quantitative map

extentextent (x1, x2, y1, y2) [μm]**features = None**Available features (see `nanite.qmap.available_features`)**get_coords**Get the qmap coordinates for each curve in *QMap.ds***Parameters which** (*str*) – “px” for pixels or “um” for microns.**group = None**Indentation data (instance of `nanite.IndentationGroup`)**shape**

shape of the map [px]

```
nanite.qmap.available_features = ['data min height', 'fit contact point', "fit young's mod"]
```

Available features for quantitative maps

List of changes in-between nanite releases.

7.1 version 1.5.5:

- setup: make tkinter optional for frozen applications

7.2 version 1.5.4:

- setup: bump scikit-learn from 0.18.0 to 0.23.0 (different model results due to bugfixes, enhancements, or random sampling procedures; the tests have been updated accordingly)
- setup: bump afmformats from 0.10.0 to 0.10.2

7.3 version 1.5.3

- setup: new builds for Python 3.8

7.4 version 1.5.2

- enh: be more verbose when tip position cannot be computed
- setup: bump afmformats from 0.7.0 to 0.10.0

7.5 version 1.5.1

- setup: bump afmformats from 0.6.0 to 0.7.0 (metadata fixes)

7.6 version 1.5.0

- feat: IndentationGroup.get_enum returns a curve from an enum value
- setup: bump afmformats from 0.5.0 to 0.6.0 (hdf5 export, improved tab export)

7.7 version 1.4.1

- enh: set parameter *baseline* to “vary” for all models
- fix: make sure that *model_key* is set before *params_initial* when fitting with kwargs (otherwise, *params_initial* might reset)

7.8 version 1.4.0

- feat: add function *Indentation.get_rating_parameters*
- feat: compute additional ancillary parameter “Maximum indentation”
- feat: new functions *model.get_parm_unit* and updated *model.get_parm_name* to work with ancillary parameters as well

7.9 version 1.3.0

- feat: allow to define ancillary parameters for models and use them during fitting by default
- feat: *Indentation.get_initial_fit_parameters* now automatically computes common and model-related ancillary parameters if no initial parameters are present
- enh: allow to set the *model_key* in more functions of *Indentation*
- ref: use *idnt* to represent Indentation instances
- fix: preprocessing steps not stored in *Indentation.preprocessing*
- setup: bump afmformats from 0.4.1 to 0.5.0

7.10 version 1.2.4

- enh: update boundaries and default values for model parameters

7.11 version 1.2.3

- fix: FitProperties did not detect changes in “params_initial”

7.12 version 1.2.2

- setup: bump afmformats version from 0.3.0 to 0.4.1

7.13 version 1.2.1

- enh: skip computation of tip position if it is already in the dataset and cannot be computed e.g. due to missing spring constant
- fix: typo in `get_data_paths_enum`
- setup: bump afmformats version from 0.2.0 to 0.3.0

7.14 version 1.2.0

- tests: `np.asscalar` is deprecated
- ref: migrate to afmformats (#1)
- docs: minor improvements

7.15 version 1.1.2

- fix: add `__version__` property
- tests: use `time.perf_counter` for timing tests
- docs: improved LaTeX rendering

7.16 version 1.1.1

- setup: migrate to PEP 517 (`pyproject.toml`)
- docs: minor update

7.17 version 1.1.0

- feat: add contact point to available features in qmap visualization
- fix: avoid two invalid operations when computing features

7.18 version 1.0.1

- fix: invalid operation when loading data with a callback function

7.19 version 1.0.0

- docs: minor update

7.20 version 0.9.3

- enh: store nanite and h5py library versions in rating container
- enh: update hyperparameters of rating regressors
- ref: deprecation in h5py: replace `dataset.value` by `dataset[...]`

7.21 version 0.9.2

- ref: renamed the mode *model_hertz_parabolic* to *model_hertz_paraboloidal* to be consistent
- docs: updat code reference and other minor improvements

7.22 version 0.9.1

- fix: *preprocessing* keyword not working in *Indentation.fit_model*
- docs: add another scripting example and minor improvements
- tests: increase coverage

7.23 version 0.9.0

- ref: remove legacy “discrete” feature type
- ref: renamed kwargs for *Indetation.rate_quality*
- ref: new method *nanite.load_group* for loading experimental data
- ref: new class `read.data.IndentationData` for managing data
- ref: replace `dataset.IndentationDataSet` with `group.IndentationGroup` to avoid ambiguities
- fix: add missing “zef18” training set
- fix: sample weight computation failed when a rating level was missing
- enh: add *nanite-generate-training-set* command line program
- tests: reduce warnings and increase coverage
- cleanup: old docs in `nanite.rate.io`
- docs: major update using helper extensions

7.24 version 0.8.0

- initial release

CHAPTER 8

Bibliography

CHAPTER 9

Indices and tables

- `genindex`
- `modindex`
- `search`

Bibliography

- [Bil92] G. G. Bilodeau. Regular pyramid punch problem. *Journal of Applied Mechanics*, 59(3):519, 1992. doi:10.1115/1.2893754.
- [Dob18] Wolfgang Dobler. Truncated power series approximation for the relationship between indentation and force of a spherical indenter in atomic force microscopy. personal communication with Wolfgang Dobler, JPK Instruments, Berlin, November 2018.
- [Lov39] A. E. H. Love. Boussinesq's problem for a rigid cone. *The Quarterly Journal of Mathematics*, os-10(1):161–175, 1939. doi:10.1093/qmath/os-10.1.161.
- [MKH+19] Stephanie Möllmert, Maria A. Kharlamova, Tobias Hoche, Anna V. Taubenberger, Shada Abuhattum, Veronika Kuscha, Michael Brand, and Jochen Guck. Zebrafish spinal cord repair is accompanied by transient tissue stiffening. (*manuscript in preparation*), 2019.
- [MAM+19] Paul Müller, Shada Abuhattum, Stephanie Möllmert, Elke Ulbricht, Anna V. Taubenberger, and Jochen Guck. Nanite: using machine learning to assess the quality of atomic force microscopy-enabled nano-indentation data. *BMC Bioinformatics*, 20(1):1–9, sep 2019. doi:10.1186/s12859-019-3010-3.
- [MMG18] Paul Müller, Stephanie Möllmert, and Jochen Guck. Atomic force microscopy indentation data of zebrafish spinal cord sections. *Figshare*, 11 2018. doi:10.6084/m9.figshare.7297202.v1.
- [Sne65] Ian N. Sneddon. The relation between load and penetration in the axisymmetric boussinesq problem for a punch of arbitrary profile. *International Journal of Engineering Science*, 3(1):47–57, may 1965. doi:10.1016/0020-7225(65)90019-4.

n

- `nanite.fit`, 35
- `nanite.group`, 26
- `nanite.indent`, 23
- `nanite.model`, 28
 - `nanite.model.model_conical_indenter`, 30
 - `nanite.model.model_hertz_paraboloidal`, 31
 - `nanite.model.model_hertz_three_sided_pyramid`, 32
 - `nanite.model.model_sneddon_spherical`, 33
 - `nanite.model.model_sneddon_spherical_approximation`, 34
- `nanite.preproc`, 27
- `nanite.qmap`, 42
- `nanite.rate.features`, 36
- `nanite.rate.io`, 40
- `nanite.rate.rater`, 39
- `nanite.rate.regressors`, 40
- `nanite.read`, 27

A

ANCILLARY_COMMON (in module *nanite.model*), 29
 append() (*nanite.group.IndentationGroup* method), 26
 apply() (*nanite.preproc.IndentationPreprocessor* static method), 27
 apply_preprocessing() (*nanite.indent.Indentation* method), 23
 available() (*nanite.preproc.IndentationPreprocessor* static method), 27
 available_features (in module *nanite.qmap*), 43
 available_preprocessors (in module *nanite.preproc*), 28

C

CannotSplitWarning, 27
 compute_emodulus_mindelta() (*nanite.indent.Indentation* method), 23
 compute_emodulus_vs_mindelta() (*nanite.fit.IndentationFitter* method), 36
 compute_features() (*nanite.rate.features.IndentationFeatures* static method), 36
 compute_opt_mindelta() (*nanite.fit.IndentationFitter* static method), 36
 compute_sample_weight() (*nanite.rate.rater.IndentationRater* static method), 39
 compute_tip_position() (*nanite.preproc.IndentationPreprocessor* static method), 27
 contact_point (*nanite.rate.features.IndentationFeatures* attribute), 38
 correct_force_offset() (*nanite.preproc.IndentationPreprocessor* static method), 27
 correct_split_approach_retract() (*nanite.preproc.IndentationPreprocessor* static method), 27

correct_tip_offset() (*nanite.preproc.IndentationPreprocessor* static method), 28

D

data (*nanite.indent.Indentation* attribute), 26
 datafit_apr (*nanite.rate.features.IndentationFeatures* attribute), 38
 DataMissingWarning, 42
 datares_apr (*nanite.rate.features.IndentationFeatures* attribute), 38
 dataset (*nanite.rate.features.IndentationFeatures* attribute), 38
 datasets (*nanite.rate.io.RateManager* attribute), 41
 datax_apr (*nanite.rate.features.IndentationFeatures* attribute), 38
 datay_apr (*nanite.rate.features.IndentationFeatures* attribute), 39
 delta_of_a() (in module *nanite.model.model_sneddon_spherical*), 33

E

estimate_contact_point_index() (*nanite.indent.Indentation* method), 24
 estimate_optimal_mindelta() (*nanite.indent.Indentation* method), 24
 export() (*nanite.indent.Indentation* method), 24
 export_training_set() (*nanite.rate.io.RateManager* method), 40
 extent (*nanite.qmap.QMap* attribute), 42

F

feat_bin_apr_spikes_count() (*nanite.rate.features.IndentationFeatures* method), 37
 feat_bin_cp_position() (*nanite.rate.features.IndentationFeatures* method), 37

- feat_bin_size() (*nanite.rate.features.IndentationFeatures* method), 37
 - feat_con_apr_flatness() (*nanite.rate.features.IndentationFeatures* method), 37
 - feat_con_apr_size() (*nanite.rate.features.IndentationFeatures* method), 37
 - feat_con_apr_sum() (*nanite.rate.features.IndentationFeatures* method), 37
 - feat_con_bln_slope() (*nanite.rate.features.IndentationFeatures* method), 37
 - feat_con_bln_variation() (*nanite.rate.features.IndentationFeatures* method), 37
 - feat_con_cp_curvature() (*nanite.rate.features.IndentationFeatures* method), 37
 - feat_con_cp_magnitude() (*nanite.rate.features.IndentationFeatures* method), 37
 - feat_con_idt_maxima_75perc() (*nanite.rate.features.IndentationFeatures* method), 38
 - feat_con_idt_monotony() (*nanite.rate.features.IndentationFeatures* method), 38
 - feat_con_idt_spike_area() (*nanite.rate.features.IndentationFeatures* method), 38
 - feat_con_idt_sum() (*nanite.rate.features.IndentationFeatures* method), 38
 - feat_con_idt_sum_75perc() (*nanite.rate.features.IndentationFeatures* method), 38
 - feat_data_min_height_measured_um() (*nanite.qmap.QMap* method), 42
 - feat_fit_contact_point() (*nanite.qmap.QMap* method), 42
 - feat_fit_youngs_modulus() (*nanite.qmap.QMap* method), 42
 - feat_meta_rating() (*nanite.qmap.QMap* method), 42
 - feat_meta_scan_order() (*nanite.qmap.QMap* method), 42
 - features (*nanite.qmap.QMap* attribute), 43
 - fit() (*nanite.fit.IndentationFitter* method), 36
 - fit_model() (*nanite.indent.Indentation* method), 24
 - fit_properties (*nanite.indent.Indentation* attribute), 26
 - FitDataError, 35
 - FitKeyError, 35
 - FitProperties (class in *nanite.fit*), 35
 - FitWarning, 35
- ## G
- get_a() (in module *nanite.model.model_sneddon_spherical*), 33
 - get_anc_parm_keys() (in module *nanite.model*), 28
 - get_anc_parms() (in module *nanite.model*), 28
 - get_ancillary_parameters() (*nanite.indent.Indentation* method), 25
 - get_available_training_sets() (in module *nanite.rate.rater*), 40
 - get_coords (*nanite.qmap.QMap* attribute), 43
 - get_cross_validation_score() (*nanite.rate.io.RateManager* method), 41
 - get_data_paths() (in module *nanite.read*), 27
 - get_data_paths_enum() (in module *nanite.read*), 27
 - get_enum() (*nanite.group.IndentationGroup* method), 26
 - get_feature_funcs() (*nanite.rate.features.IndentationFeatures* class method), 38
 - get_feature_names() (*nanite.rate.features.IndentationFeatures* class method), 38
 - get_init_parms() (in module *nanite.model*), 28
 - get_initial_fit_parameters() (*nanite.indent.Indentation* method), 25
 - get_initial_parameters() (*nanite.fit.IndentationFitter* method), 36
 - get_model_by_name() (in module *nanite.model*), 28
 - get_parm_name() (in module *nanite.model*), 29
 - get_parm_unit() (in module *nanite.model*), 29
 - get_qmap() (*nanite.qmap.QMap* method), 42
 - get_rater() (in module *nanite.rate.rater*), 40
 - get_rates() (*nanite.rate.io.RateManager* method), 41
 - get_rating_parameters() (*nanite.indent.Indentation* method), 25
 - get_training_set() (*nanite.rate.io.RateManager* method), 41
 - get_training_set_path() (*nanite.rate.rater.IndentationRater* static method), 39
 - group (*nanite.qmap.QMap* attribute), 43
 - guess_initial_parameters() (in module *nanite.fit*), 36

H

has_contact_point
(*nanite.rate.features.IndentationFeatures*
attribute), 39

hash_file (in module *nanite.rate.io*), 42

hdf5_rated() (in module *nanite.rate.io*), 41

hertz_conical() (in module
nanite.model.model_conical_indenter), 30

hertz_paraboloidal() (in module
nanite.model.model_hertz_paraboloidal),
31

hertz_sneddon_spherical_approx() (in mod-
ule *nanite.model.model_sneddon_spherical_approximation*),
34

hertz_spherical() (in module
nanite.model.model_sneddon_spherical),
33

hertz_three_sided_pyramid() (in module
nanite.model.model_hertz_three_sided_pyramid),
32

I

Indentation (class in *nanite.indent*), 23

IndentationFeatures (class in
nanite.rate.features), 36

IndentationFitter (class in *nanite.fit*), 35

IndentationGroup (class in *nanite.group*), 26

IndentationPreprocessor (class in
nanite.preproc), 27

IndentationRater (class in *nanite.rate.rater*), 39

index() (*nanite.group.IndentationGroup* method), 26

is_fitted (*nanite.rate.features.IndentationFeatures*
attribute), 39

is_valid (*nanite.rate.features.IndentationFeatures* at-
tribute), 39

L

load() (in module *nanite.rate.io*), 41

load_data() (in module *nanite.read*), 27

load_group() (in module *nanite.group*), 26

load_hdf5() (in module *nanite.rate.io*), 41

load_training_set()
(*nanite.rate.rater.IndentationRater* class
method), 39

M

meta (*nanite.rate.features.IndentationFeatures* at-
tribute), 39

model_doc (*nanite.model.nanite.model.model_submodule*
attribute), 29

model_key (*nanite.model.nanite.model.model_submodule*
attribute), 29

model_name (*nanite.model.nanite.model.model_submodule*
attribute), 29

ModelIncompleteError, 28

N

names (*nanite.rate.rater.IndentationRater* attribute), 40

nanite.fit (module), 35

nanite.group (module), 26

nanite.indent (module), 23

nanite.Indentation (built-in class), 23

nanite.IndentationGroup (built-in class), 23

nanite.IndentationRater (built-in class), 23

nanite.load_group() (built-in function), 23

nanite.model (module), 28

nanite.model.model_conical_indenter
(module), 30

nanite.model.model_hertz_paraboloidal
(module), 31

nanite.model.model_hertz_three_sided_pyramid
(module), 32

nanite.model.model_sneddon_spherical
(module), 33

nanite.model.model_sneddon_spherical_approximation
(module), 34

nanite.model.model_submodule.get_parameter_defaults
(in module *nanite.model*), 29

nanite.model.model_submodule.model() (in
module *nanite.model*), 29

nanite.model.model_submodule.residual()
(in module *nanite.model*), 29

nanite.preproc (module), 27

nanite.QMap (built-in class), 23

nanite.qmap (module), 42

nanite.rate.features (module), 36

nanite.rate.io (module), 40

nanite.rate.rater (module), 39

nanite.rate.regressors (module), 40

nanite.read (module), 27

O

obj2str() (in module *nanite.fit*), 36

P

parameter_keys (*nanite.model.nanite.model.model_submodule*
attribute), 29

parameter_names (*nanite.model.nanite.model.model_submodule*
attribute), 29

parameter_units (*nanite.model.nanite.model.model_submodule*
attribute), 30

path (*nanite.rate.io.RateManager* attribute), 41

pipeline (*nanite.rate.rater.IndentationRater* at-
tribute), 40

preprocessing (*nanite.indent.Indentation* attribute),
26

Q

QMap (*class in nanite.qmap*), 42

R

rate() (*nanite.rate.rater.IndentationRater method*), 40

rate_quality() (*nanite.indent.Indentation method*),
25

RateManager (*class in nanite.rate.io*), 40

ratings (*nanite.rate.io.RateManager attribute*), 41

reg_names (*in module nanite.rate.regressors*), 40

reg_trees (*in module nanite.rate.regressors*), 40

register_model() (*in module nanite.model*), 29

reset() (*nanite.fit.FitProperties method*), 35

reset() (*nanite.indent.Indentation method*), 26

restore() (*nanite.fit.FitProperties method*), 35

S

samples (*nanite.rate.io.RateManager attribute*), 41

save_hdf5() (*in module nanite.rate.io*), 41

shape (*nanite.qmap.QMap attribute*), 43

smooth_height() (*nanite.preproc.IndentationPreprocessor
static method*), 28

subset_with_path()
(*nanite.group.IndentationGroup method*),
26

V

VALID_FEATURE_TYPES (*in module
nanite.rate.features*), 39

verbose (*nanite.rate.io.RateManager attribute*), 41